




# LABORATORY LINES

ICE MAKER



## Kalstein®

A different accompaniment, at your service

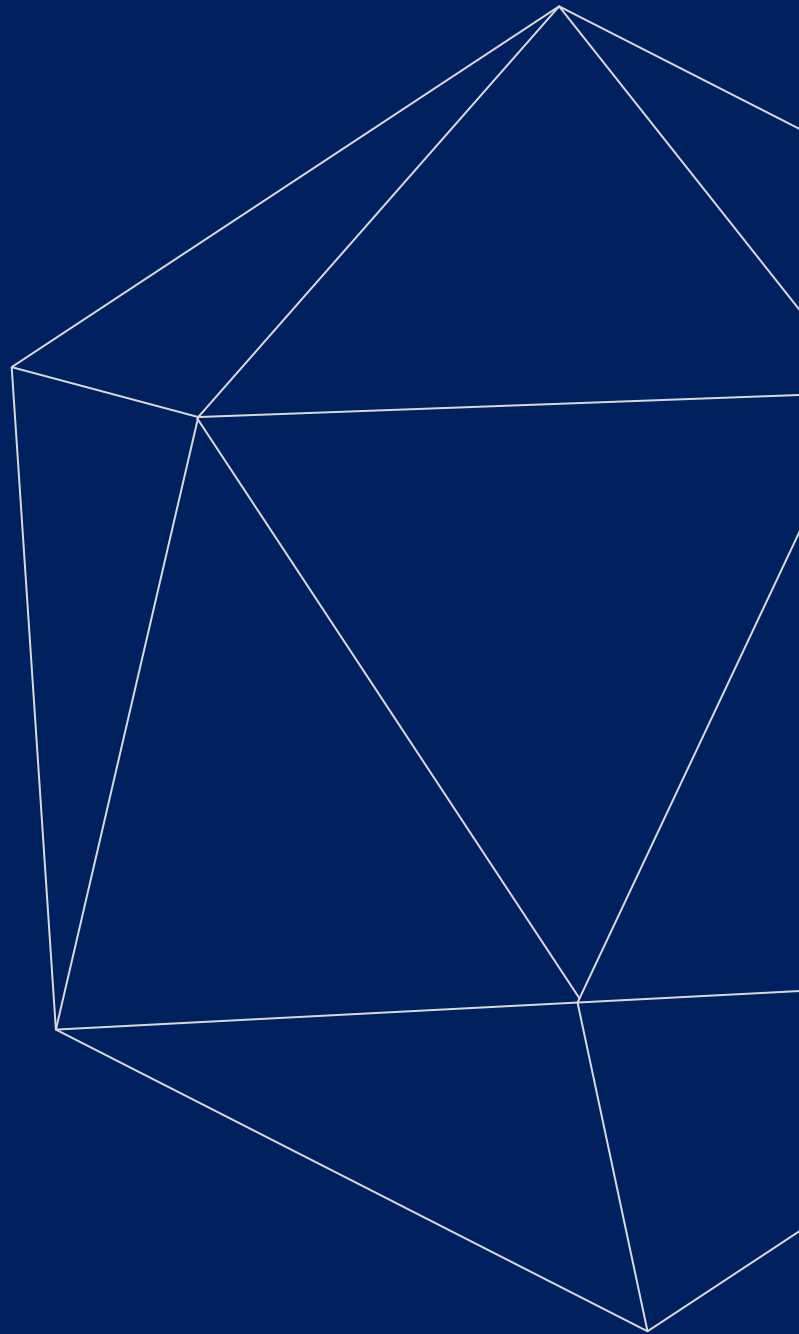



# Kalstein®

A different accompaniment, at your service

OUR SERVICES

## Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



### Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



### Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



### #Letsgivemore ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



### Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



### Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



### Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.





YR05132 / YR05133



YR05132 / YR05133

## Flake Ice Maker Machine

Serie YR05132 / YR05133

### Features:

- ✓ High-quality stainless steel is used anti-corrosion and durable. The compact structure is simple and saves space.
- ✓ The cabinet is filled with fluorine-free thermal foam with good insulation effect. The inner lining is food grade plastic. They save energy and are respectful with the environment.
- ✓ World famous brand CFC-free R134a compressors are used. The main parts and components have passed the relevant safety certifications. The product has stable quality and long service life.
- ✓ The ice making process is fully computer controlled, using imported computer chips to make the control reliable and easy to operate.
- ✓ Advanced technologies and imported accessories make low noise, smooth and reliable operation. The vents and fans are deployed to ensure that the geared motor operates reliably.
- ✓ The patented technology of the separate evaporator chip ice chamber line provides high cooling efficiency and large ice making capacity.
- ✓ Unique float type water inlet system ensures no residual water. no need for an ice cleaning process or wasted water, saving water and energy.
- ✓ Long-term performance testing and electrical safety testing will be conducted before shipment, ensuring excellent product performance.
- ✓ Fully automatic operation with full ice screen, water shortage screen, operation screen and malfunction warning screen, etc. For power-off protection. They will stop automatically in case of ice or water shortage and start automatically when the ice is not full or the water supply is restored.

### Ice Features:

- ✓ Screw extrusion hob style ice making type makes compact structure and automatic water ice separation. The optimal design of the skate blade makes the ice shape small and practical.
- ✓ Flake ice cools faster than other forms of ice, which can infiltrate a narrow space, with fast cooling speed, good ice bath effect.



### #Letsgivingmore

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
Optimum Business Center 450 Rue Baden Powell,  
• 34000 Montpellier, France.  
Tél: +33 467158849 / +33 680760710/ +33 663810023  
<https://kalstein.eu>  
KALSTEIN FRANCE, S.A. S



## Flake Ice Maker Machine

Serie YR05132 / YR05133


 Energy-saving and environmental friendly



YR05132 / YR05133



YR05132 / YR05133

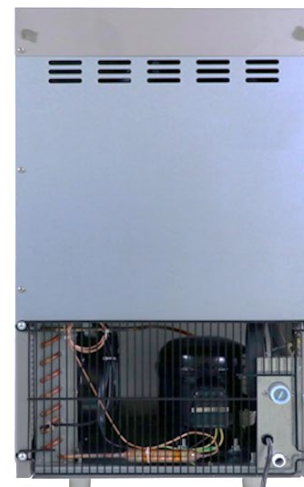

YR05132 / YR05133

### Parameters per Model

Model	YR05132	YR05133
Ice Making Capacity	40 Kg / 24 h	50 Kg / 24 h
Ice Storage	15 Kg	
Voltage	220 V/50Hz o 110 V/60 HZ (you have to pay an additional charge)	
Power	360 W	380 W
Cooling Type	Air cooling	
Refrigerant	R134a	
Water Supply	Tap water	
Net Weight (Kg)	36	
Gross Weight (Kg)	39	
Product Dimension	400 x 540 x 700mm	
Packaging Dimension	470 x 610 x 750mm	
Stuffing Qty	20'GP	135
	40'GP	285
	40'HQ	285



#Letsgivingore 

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
Optimum Business Center 450 Rue Baden Powell,  
• 34000 Montpellier, France.  
Tél: +33 467158849 / +33 680760710/ +33 663810023  
<https://kalstein.eu>  
KALSTEIN FRANCE, S.A. S





## Flake Ice Maker Machine

Serie YR05132 / YR05133

01

### Thickened Insulation Layer


Strong insulation bin, ice can be stored for 5-8 hours after ice full shutdown.



02

### SUS304

High-quality 304 stainless steel, beautiful and durable.

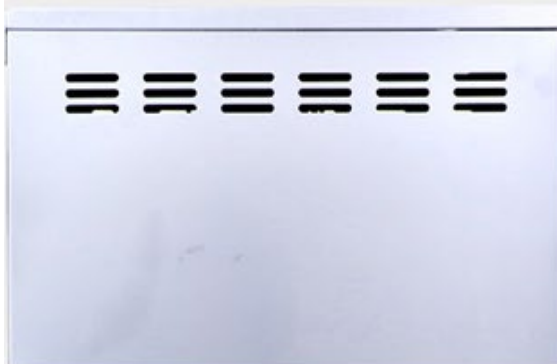

03

### Water Filter

Set the filter at the inlet, to prevent large particles into the machine.



04

### All-round Cooling

Cooling holes are installed around the machine, to avoid high temperature operation of the compressor.



#Let'sgive more ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
Optimum Business Center 450 Rue Baden Powell,  
• 34000 Montpellier, France.  
Tél: +33 467158849 / +33 680760710/ +33 663810023  
<https://kalstein.eu>  
KALSTEIN FRANCE, S.A. S





YR05132 / YR05133

## Flake Ice Maker Machine

Serie YR05132 / YR05133



YR05132 / YR05133



**#Let'sgivemore** 


Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
Optimum Business Center 450 Rue Baden Powell,  
• 34000 Montpellier, France.  
Titl: +33 467158849 / +33 680760710 / +33 663810023  
<https://kalstein.eu>  
**KALSTEIN FRANCE, S.A. S**



**Articles Kalstein**  
(Social Networks and Official Websites)

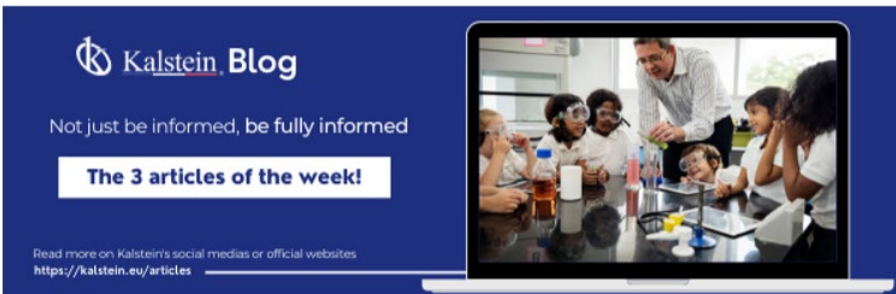

**Kalstein Blog**  
Not just be informed, be fully informed  
**The 3 articles of the week!**  
Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>



**1**  
Read this ARTICLE NOW  
**How does the climate affect your microscope?**  
A microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses and mirrors, which form...

**Kalstein Blog**  
**Article of the week:**  
How does the climate affect your microscope?  
Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>



**2**  
Read this ARTICLE NOW  
**What is the importance of the laboratory refrigerator?**  
A refrigerator is a device that is one of the most common pieces of laboratory equipment. It is used to maintain, in a controlled environment (refrigerated space), various fluids and substances...

**Kalstein Blog**  
**Article of the week:**  
What is the importance of the laboratory refrigerator?  
Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>

**Kalstein Training Videos**  
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



MONTPELLIER - FRANCE  
**HEADQUARTER**  
Optimum Business Center,  
450 Rue Baden Powell,  
34000 Montpellier - France.

**Fax:** +33 (0) 4 67 15 88 49  
**Tel:** +33 (0) 6 80 76 07 10  
[sales@kalstein.eu](mailto:sales@kalstein.eu)

<https://kalstein.eu/>  
<https://www.kalstein.fr/>

**Official Websites for Latin America.**

