


# Kalstein A different accompaniment, at your service

## LABORATORY LINES ELECTROPHORESIS



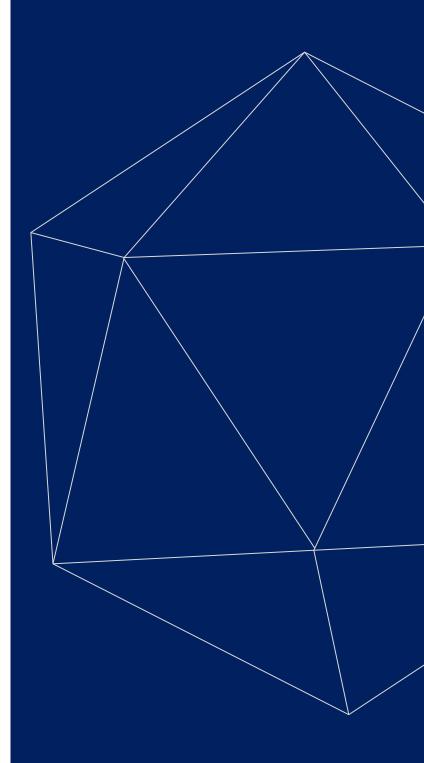

All rights reserved ® KALSTEIN France S. A. S., 2 Rue Jean Lantier • 75001 Paris • +33 1 78 95 87 89 / +33 6 80 76 07 10 • https://kalstein.eu KALSTEIN FRANCE, S. A. S





**OUR SERVICES** 

#### **Benefits and Support**

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



### Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



#### **Quick Response**

Our work team is always available to response all your consults or questions, in order to support you in any situation.





#### #Letsgivemore 💗

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



#### **Technical Support**

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.





#### **Delivery Logistics**

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or.air.



#### Kalstein Worldwide


With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



#### **Horizontal Electrophoresis Tank**

Model YR03421



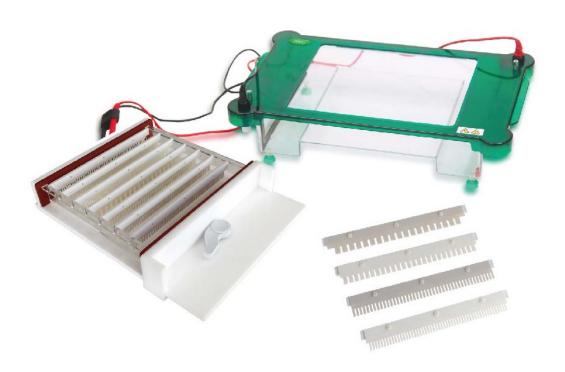

- Abundant buffer not only guarantees cooling effects, but also keeps the pH value stable during the whole experiment process.
- Removable electrodes enable maintance and replacement easy and quick.
- Adding sample more quick and suitable for the multi-channel pipettes.
- With a special gel casting, preventing the gel leakage, operation is simple and convenient.
- It's more convenient to open the lid with the easy-open buttons.
- It prevents buffer volatilizing and electricity leakage with transparent lid.
- Auto-switch-off when the lid is moved.

#### **Model YR03421 - Technical Specifications**

Gel size (W×L)	245×250 (mm)
No. of samples	19~364
Total buffer volume	2000 (ml)
Dimension (L×W×H)	487×300×93 (mm)
Weight	3.8 (kg)

#### Model YR03421 - Packing List

Parts name	Quantity
Lid with cables	1
Body tank with electrode	1
245×250mm Gel tray	1
Gel casting stand	1
1.0mm thickness 19-well comb	2
1.8mm thickness 19-well comb	2
1.0mm thickness 26-well comb	2
1.8mm thickness 26-well comb	2
1.0mm thickness 44-well comb	2
1.8mm thickness 44-well comb	2
1.0mm thickness 52-well comb	2
1.8mm thickness 52-well comb	2
level adjusting button	4



#### #Letsgivemore 💗

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved ® KALSTEIN France S. A. S. 2 Rue Jean Lantier • 75001 Paris • +33 1 78 95 87 89 /+33 6 80 76 07 10 • https://kalstein.eu

KALSTEIN FRANCE, S. A. S




KALSTEIN ELECTROPHORESIS SYSTEMS

## In Kalstein you can find the ideal Electrophoresis systems for your laboratory

Electrophoresis is a common technique in the clinical laboratory. It allows the separation of chemical species (nucleic acids or proteins) along an electric field based on their size and their electric charge.



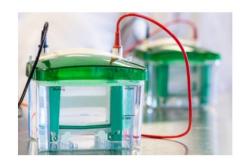

#### Agarose gel electrophoresis

Agarose gel is a good electrophoretic support. In addition, its nature gives it the characteristic of being able to separate molecules based on their size. This support allows nucleic acids to be separated based on their size. It is used to separate large molecules.



#### Polyacrylamide gel electrophoresis

Polyacrylamide gel is considered the best support. They also have great mechanical stability and are insoluble in water.

Among so many advantages there is a drawback that is limiting its use: it is neurotoxic. It is used for the separation of proteins and, like the agarose gel, allows the separation of molecules based on their size.



#### Capillary electrophoresis

This technique differs a lot from the previous ones. The support is different, since it uses fused silica capillaries, covered with a layer of polyimine. The results are visualized through a computer screen.

KALSTEIN UPDATED

## What is your ideal electrophoresis systems for your laboratory?

There are countless models, so it is normal that you do not know which electrophoresis system to buy to suit your needs. At Kalstein, we evaluate them to find what you are looking for.







#### #Letsgivemore 💗

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved ® KALSTEIN France S. A. S., 2 Rue Jean Lantier • 75001 Paris • +33 1 78 95 87 89 /+33 6 80 76 07 10 • https://kalstein.eu

KALSTEIN FRANCE, S. A. S







#LETSGIVEMORE >

## WITH THE ADCQUISITION OF A KALSTEIN EQUIPMENT

YOU MAKE A CONTRIBUTION TO:





https://onetreeplanted.org



Solidarité et coopération biomédicale

https://www.humatem.org/



https://www.maniapurefoundation.org/



#### **ANY DOUBT?** Contact us!

#### **Articles Kalstein**

(Social Networks and Official Websites)







Article of the week: How does the climate affect your microscope?





Article of the week: What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites https://kalstein.eu/articles



PARIS - FRANCE **HEADQUARTER** 

2 Rue Jean Lantier 75001 Paris - France

**Tlf**: +33 (0) 01 78 95 87 89 **Cel**: +33 (0) 6 80 76 07 10 sales@kalstein.eu

https://kalstein.eu/ https://www.kalstein.fr/

### Kalstein Training Videos (Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago

