



LABORATORY LINES

AUTOCLAVE



Kalstein®

A different accompaniment, at your service





Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ❤️

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



LED Display Portable Pressure Steam Autoclave

Series Model YR05701



Safety Protection Device

- ✓ Safe protection of water lacking
- ✓ The buzzer alarms when a cycle of operation is completed.
- ✓ Over temperature & over pressure auto-protection.
- ✓ Over current protector.
- ✓ Safety valve for automatically released against over-pressure.

Main Configuration

- ✓ Control: PID Microprocessor
- ✓ Status display: LED display
- ✓ Steam generator: Built-in steam generator, stainless steel 304
- ✓ Chamber material: Stainless steel 304
- ✓ Loading: Stainless steel drum.
- ✓ Language: English
- ✓ Voltage: 220V single phase, 50Hz

Application

YR series portable vertical pressure steam autoclave is the equipment that disinfects and sterilizes the objects with saturated steam, its structure is vertical and heating way is electric heating. It is composed of sterilization room, shell, door, steam generator, pipe system, control system, display system and power-supply system.

It is widely used in CSSD, lab, scientific research institution, pharmaceutical, food industry and is suitable for items which are high temperature-resisting and high pressure-resisting including unwrapped or wrapped instrument, dressing, rubber, liquid, plant-in instrument, porous item, hollow instrument to be sterilized.

Product Features

- ✓ High Quality SUS304 stainless steel material.
- ✓ Wing nut locking mechanism.
- ✓ Door safety interlock device.
- ✓ Electric heating system.
- ✓ Double scale indicator pressure gauge.
- ✓ LED digital display temperature, pressure and time, error code.
- ✓ PID control, temperature and time can be set freely.
- ✓ High-quality silicone seal, high temperature resistant, durable and not deformed.
- ✓ Double safety valve to ensure safe.
- ✓ High pressure faucet drainage water.
- ✓ Water shortage protection.
- ✓ Easy operate, safe and reliable.
- ✓ Automatic cut off the power if the water is not enough.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
 2 Rue Jean Lantier • 75001 Paris •
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
 KALSTEIN FRANCE, S. A. S



LED Display Portable Pressure Steam Autoclave

Series Model YR05701



Model	YR05701	YR05702	YR05703
Volume	18L	24L	30L
Chamber size	Ø300×330mm	Ø300×460mm	Ø300×520mm
Working temperature	115-129°C	115-129°C	115-129°C
Working pressure	0.145-0.165Mpa	0.145-0.165Mpa	0.145-0.165Mpa
Time scope	0-999min	0-999min	0-999min
Packing size	400×400×540mm	400×400×640mm	400×400×740mm
G.W/N.W	17kg/16kg	19kg/18kg	21kg/20kg
Power	AC220V/50Hz/2.0KW	AC220V/50Hz/2.0KW	AC220V/50Hz/2.0KW



#Let'sgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale


<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

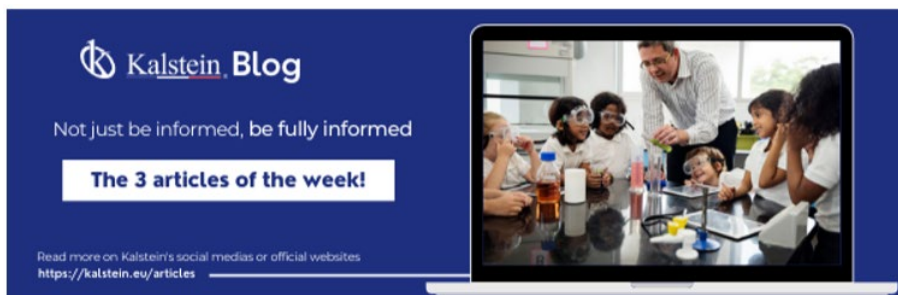

Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW

How does the climate affect your microscope?

microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses and mirrors, which form

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW

What is the importance of the laboratory refrigerator?

refrigerator is a device that is one of the most common pieces of equipment used in laboratories. It is a controlled environment (refrigerated space), various fluids and substances

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?


Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances

23 views • 3 years ago

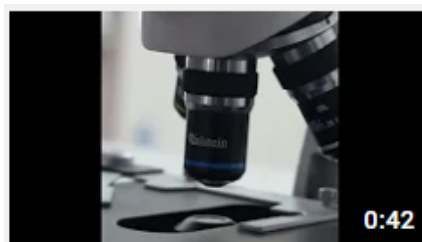

Kalstein's Microscopes

445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER

2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89

Cel: +33 (0) 6 80 76 07 10

sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

