




LABORATORY LINES

CABINETS



Kalstein®

A different accompaniment, at your service

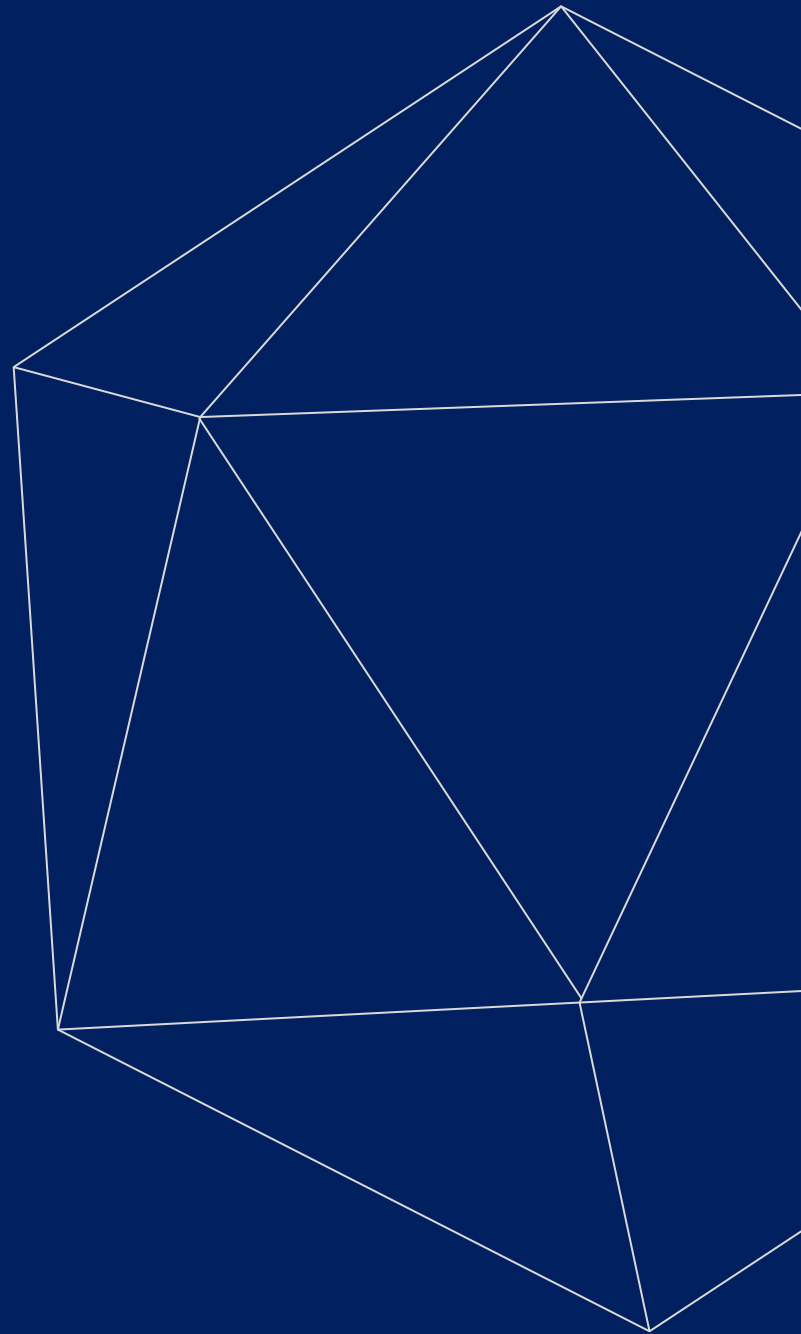



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Cabinets
Laminar Flow Cabinet
Model YR05740



- **The material of the work area is 304 stainless steel**
- **Arbitrary positioning door shifting system**
- **Digital display LIQUID crystal control interface**



Product Description

This is a kind of clean bench widely used in medical treatment, chemical experiment, electron, precise machinery, instrument, apparatus, medical preparation, operating room, foodstuff, drinks, daily use cosmetic, stock raising and aviation, as well as providing germfree dustless clean environment.

1. Adopting arbitrary positioning sliding door system
2. The shell is welded by cold-rolled plates, the surface is electrostatically sprayed, and the work surface is SUS304 brushed stainless steel, which is corrosion-resistant and easy to clean.

3. The lighting and sterilization systems are interlocked safely, and have the timing function of the ultraviolet lamp

4. Digital display LCD control interface, more user-friendly design.

5. The vertical quasi-closed tabletop and the formation of a downward flow air curtain in the operation room can effectively prevent the input of external air and the operation area is clean.

6. Equipped with a HEPA filter and a primary filter for preliminary filtration, which can effectively extend the service life of the high-efficiency filter.



#Letsgivemore ♥


Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
 Optimum Business Center 450 Rue Baden Powell,
 • 34000 Montpellier, France.
 Tit: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



KALSTEIN® Cabinets
Laminar Flow Cabinet
Model YR05740

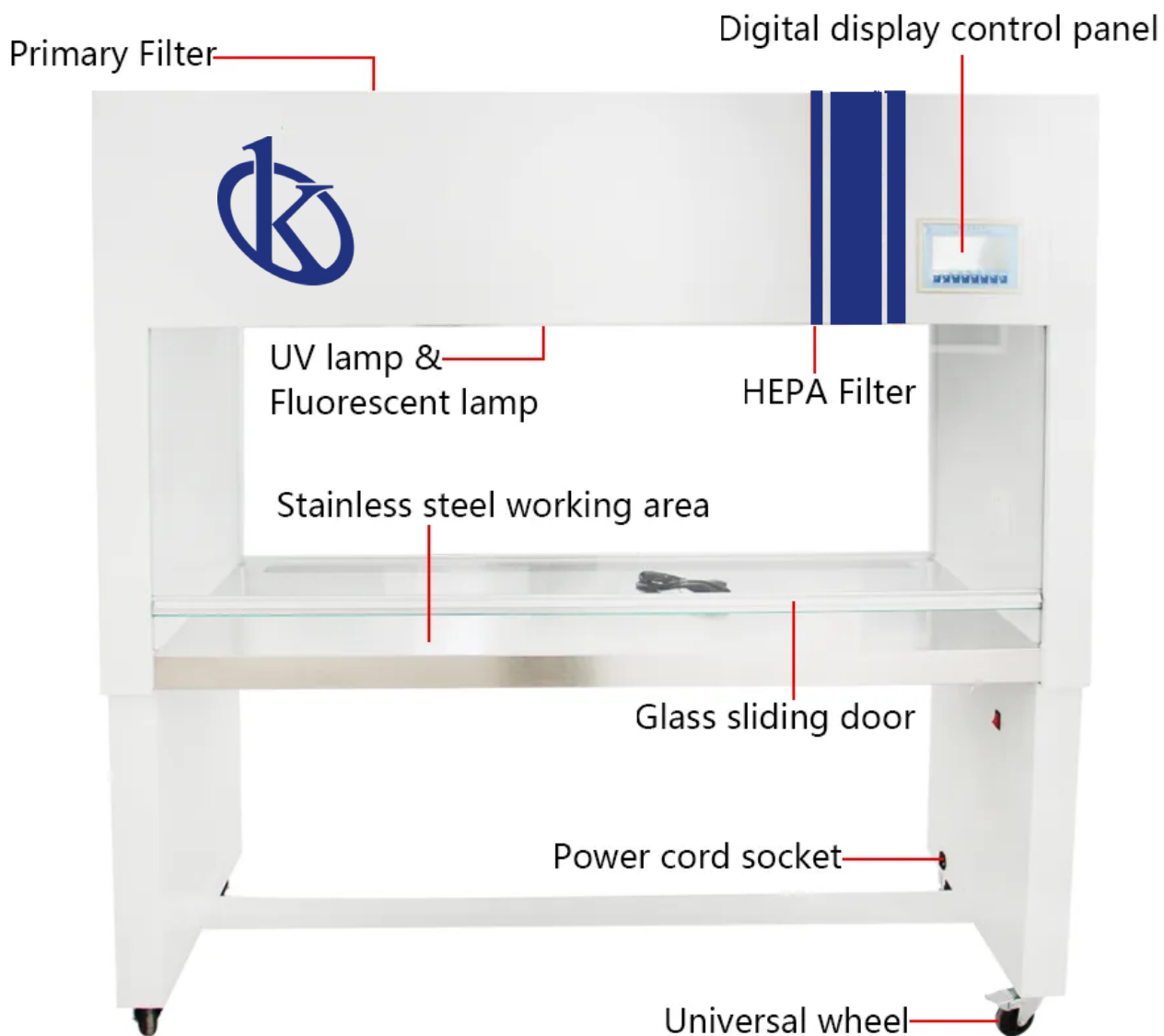

#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



KALSTEIN® Cabinets Laminar Flow Cabinet Model YR05740



Product details



High efficiency filter

HEPA Filter ensures a clean environment in the work area.

Fluorescent lamp & UV lamp

Equipped with high-brightness fluorescent lamps, safety interlocking of lighting and sterilization systems.




Socket (Four options)

The working area is equipped with sockets (four options are available) to facilitate the experiment plug-in needs.

All steel material

The product structure is firmer and the appearance is more beautiful.

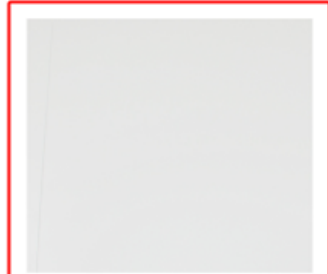


LCD control panel

Digital display control panel, three-speed adjustable for fast, medium and slow, more humanized design.

Copper core motor

The motor has strong power, low heat generation, low noise, energy-saving and power-saving.

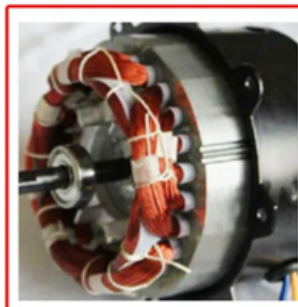

Working area

The material of the work area is 304 stainless steel, resistant to corrosion and easy to clean.

Universal wheel

The base is equipped with universal wheels to facilitate the movement of the clean bench.



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



KALSTEIN® Cabinets Laminar Flow Cabinet Model YR05740



Technical Specifications

Model	YR05740
Air supply	Vertical Air Supply
Cleanliness grade	Class 100 (US Federal 209E)
Number of colonies	≤0.5pcs/dish-hour (Ø90mm culture dish)
The average wind speed	0.25-0.45m/s
Noise	≤62dB (A)
Power supply	AC single phase 220V /50Hz
Vibration half peak	≤0.5μM (x, y, z direction)
Illumination	≥300LX
Max. Power consumption	800 W
Dimension of working area (W*D*H)(mm)	1370x650x580
Overall dimension (W*D*H)(mm)	1530x690x1600
Specification and quantity of high efficiency filter	1355x555x50 x 1
Specification and quantity of fluorescent lamp/UV lamp	12 W x 1 / 30 W x 1
Applicable persons	Double Person – Single Side



#Letsgivingmore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tit: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 


**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:

 **Fundación
Jacinto Convit**
<https://www.jacintoconvit.org/>


 **ONE
TREE
PLANTED**

<https://onetreepanted.org>

 **Humatem**
Solidarité et coopération
biomédicale
<https://www.humatem.org/>


**MANIAPURE
FOUNDATION USA**
<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

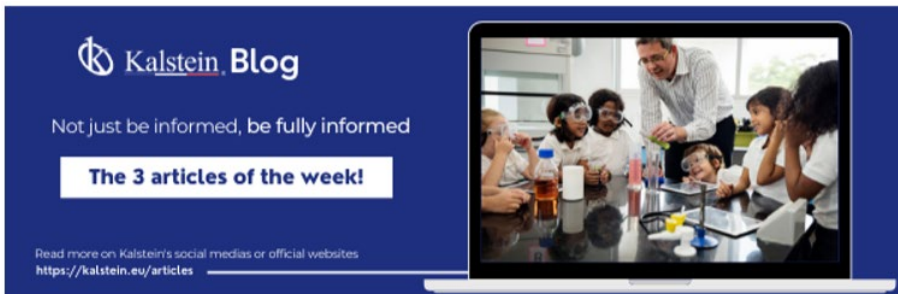

Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



PARIS - FRANCE
HEADQUARTER
2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>



1
Read this ARTICLE NOW

How does the climate affect your microscope?

microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses and mirrors, which form

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW

What is the importance of the laboratory refrigerator?

refrigerator is a device that is one of the most common pieces of equipment used in laboratories. It is a controlled environment (refrigerated space), various fluids and substances

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?


Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago

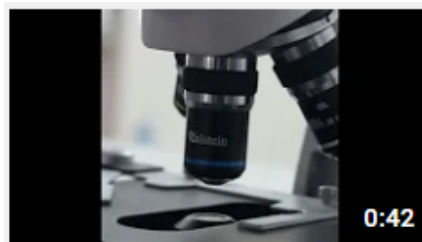

Kalstein's Microscopes

445 views • 1 week ago

