



LABORATORY SECTOR

AIR COMPRESSOR



Kalstein®

A different accompaniment, at your service





Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Air Compressor Oil Free Air Compressor Model YR0AC-30



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Air Compressor
Oil Free Air Compressor
Model YR0AC-30



Product Description

- ✓ High quality, big flow oil free swing piston compressor as a source of power, can provide stable oil less air and avoid the oil stain into human body.
- ✓ The parts are selecting high quality components and the machine made acc.
- ✓ To the requirement of clinic, focus on the features of large flow, low noise, clean air ,smooth operation, long time continuous work, long life, full automatic control .
- ✓ In the operation process, when air tank pressure reach to the minimum or maximum value, the compressor will be automatically started or shut down.

Model	Motor		Air Displacement	Discharge Pressure	Tank		Weight	Dimensions	Noise
	W	HP			L	Gal			
YR0AC-30	600x1	0.75	45	0.75	30	6.5	30	43x43x67	56-65



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
 2 Rue Jean Lantier • 75001 Paris •
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Air Compressor Oil Free Air Compressor Model YR0AC-30



Compressor motor build in pure copper coil, not other cheap compressor's aluminum coil, long time usage no maintained..

2pcs build in high-power cooling fan for better heat radiation.

Self-designed power box, put all power line and capacitor inside, more space saving. Power switch, start switch separately design, improve the use efficiency.

Safety valve, When the pressure of the tank is greater than 8.2 kilograms, the safety valve would auto open and deflated guarantee the safety of use.

Solenoid Valve, leak extra air in the pipe and valve. Make compressor restart more smoothly.

Inside of air tank have been antibacterial coated to prevent corrosion, and promise the air delivery up to medical grade.

High pressure pipe tightly connect to the cylinder, ensure continuous operation of the air compressor.

High quality check valve, stop the air flow back.

Filter Regulator, Regulate the air out pressure and filter out air.

pure cooper valve, not cheap iron or aluminum type.



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Air Compressor Oil Free Air Compressor Model YR0AC-30



Advantage:

SUPER SILENT

- ✓ Low working noise, can be used for all medical equipment, which need quiet air supply.

LOW VIBRATION

- ✓ With special rubber feet, reduce vibration during operation.

PURE AIR FLOW

- ✓ Oil free design, no lubrication oil needed during operation.

CORE TECHNOLOGY

- ✓ Diamond hardness cylinder ensure durable working performance.

FASHION AND DURABLE DESIGN

- ✓ Compact structure, light weight; under normal situation, the compressor can be used for more than 20000hours.

USE SAFETY

- ✓ With multiple self protection system, if here will be abnormal with pressure, current or voltage, the motor would cut off automatically to ensure equipment and personal safety.

EASY OPERATION

- ✓ Quite simple operation, connect to power supply, then no need any more maintenance, just drainage regularly.

LOW ENERGY CONSUMPTION

- ✓ Full automatic design, automatic stop and restart control, low consumption

HIGH PRECISION FILTRATION

- ✓ With double filters, ensure high precision of outlet air flow

TANK INSIDE HAS DONE ANTI-RUST TREATMENT

- ✓ All tanks inside with anti-rust treatment, ensure pure outlet air flow for medical equipment



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale


<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

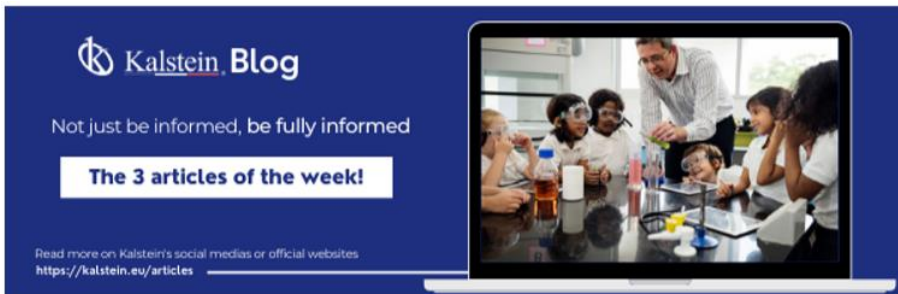

Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

1
Read this ARTICLE NOW



How does the climate affect your microscope?

microscope is an instrument that allows observing and measuring small objects. This is achieved through a series of lenses that focus light on the object being observed.

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

2
Read this ARTICLE NOW



What is the importance of the laboratory refrigerator?

Refrigerator is a device that is one of the most essential equipment that is used to maintain, in a controlled environment (refrigerated space), various fluids and substances.

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER

2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89

Cel: +33 (0) 6 80 76 07 10

sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

